nenò[®]

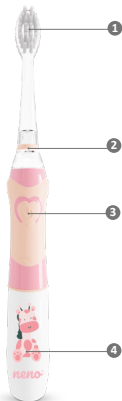
Fratelli

Table of contents

Spis treści

4	Instrukcja obsługi	PL
7	User manual	EN
10	Bedienungsanleitung	DE
14	Manuel de l'utilisateur	FR
18	Návod k použití	CZ
21	Používateľská príručka	SK
24	Használati utasítás	HU
28	Användarhandbok	SE
31	Käyttäjän käsikirja	FI
34	Brukerhåndbok	NO
38	Brugermanual	DK
41	Gebruikershandleiding	NL
45	Manual del usuario	ES
48	Manuale utente	IT
52	Manual de utilizare	RO
55	Korisnički priručnik	HR/BA
58	Uputstvo za upotrebu	RS/ME
62	Uporabniški priročnik	SI
65	Εγχειρίδιο χρήστη	GR
69	Lietotāja rokasgrāmata	LV
72	Instrukcija	LT
75	Kasutusjuhend	ET
79	Керівництво користувача	UA
82	Ръководство за употреба	BG
86	Упатство за употреба	MK
90	Manual do utilizador	PT

A**B****C**

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EN

USER MANUAL

Dear Customer, thank you for purchasing NENO FRATELLI. The product you have obtained is an electric toothbrush for small children. This device is used to brush the teeth of a baby and maintain their oral hygiene. Before using the device, please read instructions below.

01. PRECAUTIONS

1. If you do not understand the instructions for use or precautions listed in this manual, do not use the device. Incorrect use of the device may damage it.
2. Do not disassemble the product. This may result in damage to the product, electric shock and voiding the warranty.
3. Do not use accessories other than those recommended by the manufacturer.
4. Do not use flammable or explosive substances to clean the product.
5. This product is only used for brushing teeth and massaging gums. Do not use the device for any other purpose.
6. Stop using the device if it is clearly damaged.

7. Do not submerge the product in water. This can lead to a malfunction and may result in fire, electric shock or explosion.
8. Keep the device away from water containers such as a bath or sink.
9. Do not replace the battery of the device when your hands are wet.
10. Do not use used, old batteries in your device.
11. Before each use, make sure the toothbrush is securely tightened in the lower battery cover and there is a seal between the toothbrush and the battery cover.
12. Pay attention to the poles of the battery when placing it in the device.
13. Do not place the device in rooms with high moisture or high temperatures, such as bathrooms.
14. Remove the battery from the device if you will not use it for an extended period of time to avoid the risk of fire.
15. If the toothpaste used contains peroxides or bicarbonate (often present in tooth whitening pastes), wash the device thoroughly after use to avoid destroying its plastic casing.
16. If using the device causes your child discomfort or bleeding from the gums, consult your dentist before continuing to use it.
17. Children may only use the appliance under supervision adults who have read the operating instructions.
18. The device is not a toy – keep it away from children, so that they do not use it without the supervision of an adult.
19. The device is used for personal hygiene and should not be used on many people in dental offices.

02. DESCRIPTION OF THE DEVICE

SEE FIGURE A

1. Head with bristles
2. Light
3. Power switch
4. Screw cap/ battery space
5. Suction cup

03. PACKAGE CONTENTS

1. Toothbrush
2. Additional interchangeable head
3. User Manual

04. POWER

The device needs 1 pieces of AAA battery for proper operation. Unscrew the lid at the bottom of the device to access the battery space. Put in the AAA battery as shown in the **FIGURE B**.

Place the lid back in place and turn it, making sure it stays in place. The bump on the lid should be under the bottom of the arrow on the housing before turning, and under the dot after turning the lid. **SEE FIGURE C**

05. USING THE DEVICE

1. Place the head with the toothbrush on the plastic part of the holder with the light, then turn the head to the left to lock it in place.
2. Press the power button of the brush to turn on the light of the device head.
3. Moisten the bristles at the toothbrush head.
4. Squeeze the toothpaste onto the toothbrush head.
5. Press the power button to switch on the brush vibration and start cleaning your teeth. The first time you use it, select 1 of the 4 modes below. The toothbrush remembers the last mode and you don't need to select it again when you start up again. When you press the power button while brushing, the toothbrush will start a new 2 minute cycle for the next mode.

Mode 1: White light and vibration (Daily use)

Mode 2: Coloured light and vibration (Make brushing a pleasure)

Mode 3: Vibration (Economical mode)

Mode 4: White light (Watch your teeth)

6. Thoroughly clean every part of your child's teeth. Every 30 seconds of brushing, the toothbrush will stop for 1 second to allow the user to change the brushed area. Correct tooth brushing should take a total of 2 minutes. The device will automatically switch off after 2 minutes.
7. When you're done brushing your teeth, turn off the toothbrush by pressing the power button again.
8. After cleaning the baby's teeth, wash the toothbrush under running water. Remove the remaining water from the bristles of the brush and use soft cloth to wipe the handle of the device.
9. When you are not using the brush, keep it upright.

NOTE: To avoid toothpaste splashing around, turn off the device before removing it from the baby's mouth. For the sake of your child's hygiene and safety, replace the brush head every three months.

06. CLEANING THE BRUSH

1. Wash the brush head and its bristles after using the device and let it dry naturally.
2. Do not press the bristles or pull them out with your fingers. This can deform the toothbrush and destroy it.
3. Wipe the handle with a soft piece of cloth and let it dry naturally.
4. Do not submerge the handle underwater.

5. Do not use water above 40 °C or acidic cleaning agents to clean the device.

07. FAQ

Q: Vibrations weaken during use:

A: You push the toothbrush too hard. It should gently touch the teeth.

08. SPECIFICATION

Power source: 1x AAA battery

Waterproof: IPX7

Nylon bristles: Yes

Strokes per minute: 24000

Built-in timer: Yes

Built-in LED light: Yes

Dimensions: 188x19x19mm

Weight: 27g

09. WARRANTY CARD

The product is covered by a 24-month warranty. The terms of the guarantee can be found at: <https://neno.pl/gwarancja>

Details, contact and website address can be found at: <https://neno.pl/kontakt>

Specifications and contents of the kit are subject to change without notice. We apologise for any inconvenience.

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EN

The crossed out trash can symbol indicates that unusable electrical or electronic devices, its accessories (such as power supplies, cords) or components (for example batteries, if included) cannot be disposed of alongside with household waste. In order to dispose of the devices or its components (for example, batteries) deliver the device to the collection point, where it will be accepted free of charge. Disposal is subject to the recast version of the WEEE Directive (2012/19/ EU) and the Directive on batteries and accumulators (2006/66 / EC). Proper disposal of the device prevents degradation of the natural environment. Information about the collection points of the facilities is issued by the competent local authorities. Incorrect disposal of waste is subject to penalties provided for by the law in force in the given area.

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