neno[®] beauty

Algo


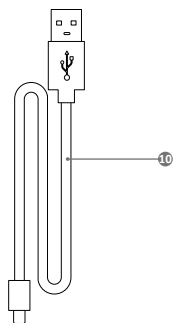
Table of contents

Spis treści

4	Instrukcja obsługi	PL
6	User manual	EN
8	Bedienungsanleitung	DE
10	Manuel de l'utilisateur	FR
12	Návod k použití	CZ
15	Používateľská príručka	SK
17	Használati utasítás	HU
19	Användarhandbok	SE
21	Käyttäjän käsikirja	FI
23	Brukerhåndbok	NO
25	Brugermanual	DK
27	Gebruikershandleiding	NL
30	Manual del usuario	ES
32	Manuale utente	IT
34	Manual de utilizare	RO
36	Korisnički priručnik	HR/BA
38	Uputstvo za upotrebu	RS/ME
40	Uporabniški pravilnik	SI
42	Εγχειρίδιο χρήστη	GR
45	Lietotāja rokasgrāmata	LV
47	Vartotojo vadovas	LT
49	Kasutusjuhend	ET
51	Посібник користувача	UA
53	Ръководство за потребителя	BG
56	Упатство за употреба	MK
58	Manual do utilizador	PT

A

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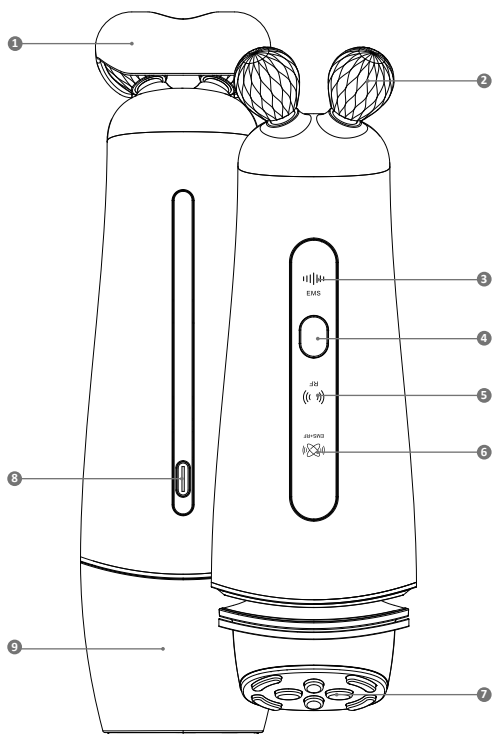Szanowny Kliencie,

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USER MANUAL

Dear Customer,

Thank you for purchasing the Neno Beauty Algo facial firming and contouring device with EMS and RF technology. The device is designed for advanced skin care - it combines EMS and radiofrequency (RF) technology with therapeutic LED light and sonic vibration to effectively firm skin, improve facial oval, stimulate microcirculation and reduce signs of aging. This user manual contains all the information necessary for proper use. Please familiarize yourself with it and keep it in case you need to use it again.

01. PRECAUTIONS

1. Do not use the device on swollen or inflamed skin, or on damaged skin, burned skin, or immediately after recent surgery.
2. The device is not recommended for people with electronic implants (e.g., pacemaker, insulin pump).
3. Do not use the device if you have metal parts in your mouth (e.g. implants, crowns, braces).
4. The device should not be used if you have active skin diseases (ulcers, psoriasis, eczema, rashes, recurrent herpes), epidermal damage (scars, wounds), hypersensitivity to electromagnetism, pregnancy and lactation, hemophilia/hemorrhagic disorders, epilepsy, atopic dermatitis (AD) and cancer.
5. Do not use the device while taking photosensitizing drugs (e.g., certain antibiotics, retinoids).
6. The device should not be used by children. Adults with physical, sensory or cognitive limitations should use the device only under the supervision of a caregiver.
7. If you have any doubts about the applicability of the device for a specific medical condition (e.g., skin diseases, any chronic illnesses, past cosmetic surgery), consult your doctor.
8. Do not apply too long on one area of the skin.
9. Avoid looking directly at operating LEDs.
10. It is recommended to remove any metal jewelry, such as earrings or chains, before using the device.
11. Avoid using massage heads around the eyes and on the neck in the thyroid region.
12. Do not use the device longer or more often than recommended. Excessive use may lead to skin irritation.
13. Do not use the device near water or in high humidity conditions, such as a bathtub, shower or sauna.
14. If you notice any adverse skin reactions during use, such as redness, itching, burning, or other unusual reactions - discontinue use immediately.
15. Transient perception of flashing lights during use is a natural phenomenon resulting from the effects of the pulses. However, if it persists afterwards, consult your doctor.
16. A slight tingling sensation felt when using the device is a normal effect of EMS pulses. If there is discomfort, it is recommended to reduce the intensity level.
17. Do not use the device while it is charging. Disconnect the device from power when it is fully charged.
18. Keep the power source away from sources of heat and water while charging.
19. Do not use the product if any component is damaged or missing.
20. Do not disassemble or repair the product yourself.
21. Do not use accessories and spare parts other than those included in the kit.
22. This product is not a toy. Do not allow children to play with the device.
23. Do not store the device and its packaging in a place easily accessible to children.
24. Do not store the device near a heat source and do not expose it directly to sunlight. Do not immerse the device in water.

02. CONTENTS OF THE SET

1. Massager

2. Storage bag
3. Microfiber cloth
4. Type-C cable
5. User Manual

03. DESCRIPTION OF THE DEVICE

SEE FIGURE A

1. Transparent protective cover
2. Diamond massage heads
3. EMS indicator
4. Button to turn on/off/change modes and intensity
5. RF indicator
6. EMS+RF indicator
7. Flat head with eight massage points
8. USB charging port
9. Magnetic cover
10. Type-C cable

04. OPERATION OF THE DEVICE

1. The device has three modes:
 - EMS: dedicated to lifting and improving skin tone, the mode uses EMS technology and diamond massage heads to precisely stimulate facial muscles, helping to firm them and increase skin elasticity. Gentle electrical pulses combined with mechanical massage enhance the facial oval and reduce signs of fatigue. The mode offers two levels of intensity, tailored to individual skin needs.
 - RF: dedicated to skin brightening and revitalization, this mode uses RF radio waves, yellow LED light (590 nm), sonic vibration and a flat head with eight massage points. It promotes blood circulation, evens out skin tone, reduces hyperpigmentation and promotes skin detoxification. RF heat stimulates cell regeneration, while sonic vibration improves the absorption of active ingredients. The mode offers two levels of intensity, allowing you to adjust to your skin's needs.
 - EMS+RF: intensive regeneration and firming mode, combining EMS technology, RF waves, red LED light (633 nm), sonic vibration and the action of a flat head with eight massage points. It strengthens facial muscles, stimulates collagen and elastin production, improves skin tone and smooths wrinkles. Red light promotes cell renewal, while sonic vibration improves circulation and absorption of active ingredients. The mode offers two levels of intensity for individual adjustment.
2. Wash and dry your face thoroughly before use. For optimal results and comfort, we recommend using the device with a conductive gel suitable for RF and EMS technology or with a light, oil-free facial product that does not impede the smooth gliding of the device over the skin.
3. To turn on the device, press and hold the A4 button (a beep will sound). Then use the same button to set the mode and intensity - the device will emit a single beep with each change. The device works in pulsed mode - a deliberate and safe action that allows for effective skin stimulation.
4. During the massage, move the device slowly across the face, making upward and outward movements - in accordance with the natural direction of lymphatic drainage. The device signals the change of the side of the face every 3 minutes, ensuring symmetrical action (beeping twice).
5. The device has a built-in 6-minute timer: the device will turn off six minutes after setting one of the modes. The time does not add up when changing intensity levels during their operation.
6. To turn off the device before 6 minutes, press and hold the A4 button - a beep will sound.
7. Once the massage is over, you can move on to the next stages of care.
8. The device has a memory mode - when the device is turned on again, it will automatically start in the last used mode.
9. Use the massager 2-4 times a week (for 6 minutes), choosing the mode and intensity level that suits your needs.

05. CLEANING AND STORAGE

1. After each use, wipe the massaging tips with a dry or slightly damp cloth.
2. Leave it to dry in a dry place, away from direct sunlight and heat sources.
3. Do not immerse the device in water and do not use highly corrosive detergents or alcohol.
4. When the product is not in use, it should be stored in the included pouch, with the magnetic lid and transparent protective cover applied.

06. DEVICE CHARGING

1. During charging, the mode indicators flash sequentially. The device cannot be used during charging.
2. The indicators of all modes light up simultaneously when the device is fully charged.
3. When the battery level is low, the running mode indicator flashes and the device beeps, indicating it needs to be recharged.
4. If you try to start the device with a fully discharged battery, it will beep four times, the mode indicators will flash, and the device will automatically turn off.

07. SPECIFICATION

Power supply: 5V 0.75A

Power: 12W

Battery: lithium-polymer 1200mAh

Timer: 6min

Material: ABS + PC + PP + zinc alloy

Dimensions: 54.5x45.3x 190 mm

Weight: 186g

08. WARRANTY CARD

The product comes with a 24-month warranty. Warranty terms and conditions can be found at: <https://neno.pl/gwarancja>

Details, contact and service address can be found at: <https://neno.pl/kontakt>

Specifications and contents are subject to change without notice. We apologize for any inconvenience.

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01. PRECAUTIONS

1. Do not use the device on swollen or inflamed skin, or on damaged skin, burned skin, or immediately after recent surgery.
2. The device is not recommended for people with electronic implants (e.g., pacemaker, insulin pump).
3. Do not use the device if you have metal parts in your mouth (e.g. implants, crowns, braces).
4. The device should not be used if you have active skin diseases (ulcers, psoriasis, eczema, rashes, recurrent herpes), epidermal damage (scars, wounds), hypersensitivity to electromagnetism, pregnancy and lactation, hemophilia/hemorrhagic disorders, epilepsy, atopic dermatitis (AD) and cancer.
5. Do not use the device while taking photosensitizing drugs (e.g., certain antibiotics, retinoids).
6. The device should not be used by children. Adults with physical, sensory or cognitive limitations should use the device only under the supervision of a caregiver.
7. If you have any doubts about the applicability of the device for a specific medical condition (e.g., skin diseases, any chronic illnesses, past cosmetic surgery), consult your doctor.
8. Do not apply too long on one area of the skin.
9. Avoid looking directly at operating LEDs.
10. It is recommended to remove any metal jewelry, such as earrings or chains, before using the device.
11. Avoid using massage heads around the eyes and on the neck in the thyroid region.
12. Do not use the device longer or more often than recommended. Excessive use may lead to skin irritation.

13. Do not use the device near water or in high humidity conditions, such as a bathtub, shower or sauna.
14. If you notice any adverse skin reactions during use, such as redness, itching, burning, or other unusual reactions - discontinue use immediately.
15. Transient perception of flashing lights during use is a natural phenomenon resulting from the effects of the pulses. However, if it persists afterwards, consult your doctor.
16. A slight tingling sensation felt when using the device is a normal effect of EMS pulses. If there is discomfort, it is recommended to reduce the intensity level.
17. Do not use the device while it is charging. Disconnect the device from power when it is fully charged.
18. Keep the power source away from sources of heat and water while charging.
19. Do not use the product if any component is damaged or missing.
20. Do not disassemble or repair the product yourself.
21. Do not use accessories and spare parts other than those included in the kit.
22. This product is not a toy. Do not allow children to play with the device.
23. Do not store the device and its packaging in a place easily accessible to children.
24. Do not store the device near a heat source and do not expose it directly to sunlight. Do not immerse the device in water.

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04. OPERATION OF THE DEVICE

1. The device has three modes:
 - EMS: dedicated to lifting and improving skin tone, the mode uses EMS technology and diamond massage heads to precisely stimulate facial muscles, helping to firm them and increase skin elasticity. Gentle electrical pulses combined with mechanical massage enhance the facial oval and reduce signs of fatigue. The mode offers two levels of intensity, tailored to individual skin needs.
 - RF: dedicated to skin brightening and revitalization, this mode uses RF radio waves, yellow LED light (590 nm), sonic vibration and a flat head with eight massage points. It promotes blood circulation, evens out skin tone, reduces hyperpigmentation and promotes skin detoxification. RF heat stimulates cell regeneration, while sonic vibration improves the absorption of active ingredients. The mode offers two levels of intensity, allowing you to adjust to your skin's needs.
 - EMS+RF: intensive regeneration and firming mode, combining EMS technology, RF waves, red LED light (633 nm), sonic vibration and the action of a flat head with eight massage points. It strengthens facial muscles, stimulates collagen and elastin production, improves skin tone and smooths wrinkles. Red light promotes cell renewal, while sonic vibration improves circulation and absorption of active ingredients. The mode offers two levels of intensity for individual adjustment.
2. Wash and dry your face thoroughly before use. For optimal results and comfort, we recommend using the device with a conductive gel suitable for RF and EMS technology or with a light, oil-free facial product that does not impede the smooth gliding of the device over the skin.
3. To turn on the device, press and hold the A4 button (a beep will sound). Then use the same button to set the mode and intensity - the device will emit a single beep with each change. The device works in pulsed mode - a deliberate and safe action that allows for effective skin stimulation.
4. During the massage, move the device slowly across the face, making upward and outward movements - in accordance with the natural direction of lymphatic drainage. The device signals the change of the side of the face every 3 minutes, ensuring symmetrical action (beeping twice).
5. The device has a built-in 6-minute timer: the device will turn off six minutes after setting one of the modes. The time does not add up when

changing intensity levels during their operation.

6. To turn off the device before 6 minutes, press and hold the A4 button - a beep will sound.
7. Once the massage is over, you can move on to the next stages of care.
8. The device has a memory mode - when the device is turned on again, it will automatically start in the last used mode.
9. Use the massager 2-4 times a week (for 6 minutes), choosing the mode and intensity level that suits your needs.

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