

neno[®]
beauty

ALVA
LED MASK
FOR FACIAL SKIN
THERAPY

Exclusive Series



blue light/ red light

Exclusive Series

Discover the Neno Beauty series

The Neno Beauty facial care series was created to accompany you in your daily rituals. Without leaving home, you gain professional care, and your skin becomes more radiant every day. With Neno Beauty, you discover your natural beauty.



Exclusive Series




Brightening mode



Anti-aging mode




Soothing mode

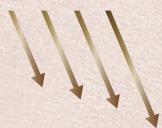


4 light wavelengths



Cordless

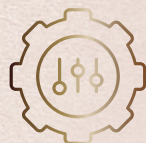

3 programs

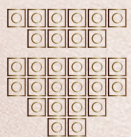

220 diodes



Timer



Exclusive Series

Four light wavelengths

470 nm

BLUE LIGHT

Blue light with a wavelength of 470 nm stimulates the skin's natural repair processes, supporting its regeneration. It has antibacterial properties, which help reduce inflammation, soothe redness and irritation, and additionally assist in regulating sebum production.

630 nm

RED LIGHT

Near-infrared light with a wavelength of 850 nm penetrates deep into the skin, stimulates tissue renewal processes, and improves firmness and smoothness, giving the complexion a fresh, youthful appearance.

660 nm

DEEP RED LIGHT

Deep red light with a wavelength of 660 nm supports collagen production, enhances skin elasticity, smooths fine lines, firms the skin, and gives the complexion a healthy glow.

850 nm


NEAR-INFRARED LIGHT

Red light with a wavelength of 630 nm stimulates microcirculation, helps even out skin tone and reduce discoloration, and also has a soothing and calming effect on the complexion.

Exclusive Series

Meet Neno Beauty Alva

Discover the power of modern skincare with the Neno Beauty Alva LED mask. Red light stimulates regeneration and improves skin elasticity, while blue light has antibacterial properties and helps fight imperfections. It is a safe and convenient therapy you can perform in the comfort of your own home.

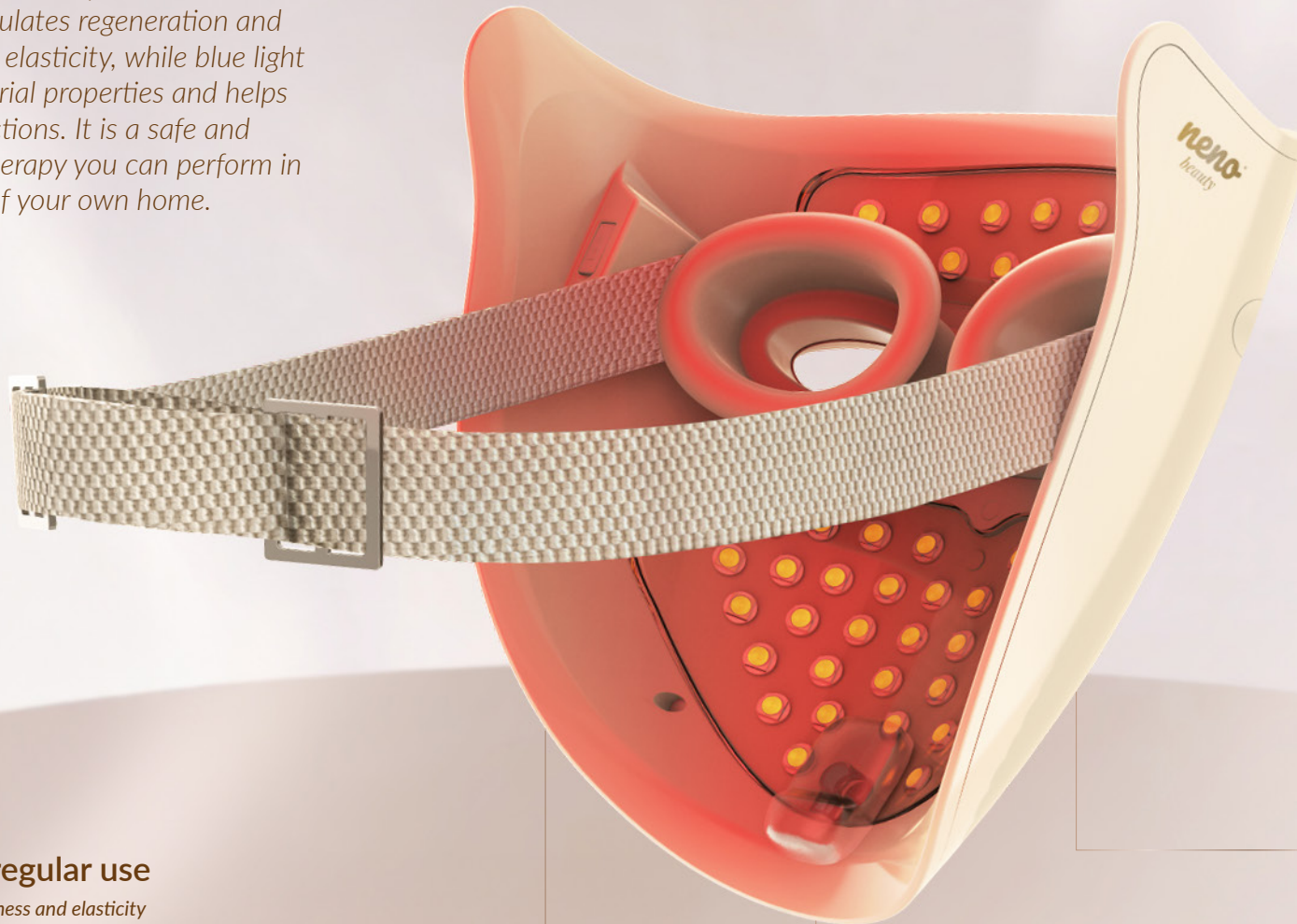

Effects of regular use

- Improved skin firmness and elasticity
- Reduction of wrinkles and fine lines
- Brightening and evening out of skin tone
- Support for skin regeneration

3 light therapy programmes

- Soothing mode
- Brightening mode
- Anti-aging mode

220 LED
Full face

Automatic timer
10 minutes

**Ultra-lightweight
and comfortable design**

High-quality silicone around the eyes,
elastic headband and chin support

Store bag

Cordless

4 precise light wavelengths

- Blue light (470 nm)
- Red light (630 nm)
- Deep red light (660 nm)
- Near infrared light (850 nm)

Easy to use

- One function button
- Voice notifications for every mode

Exclusive Series



Light therapy at home

With the Neno Beauty Alva mask, you can enjoy light therapy without leaving your house. At-home phototherapy means convenience, time-saving, and comfort, while also allowing you to care for your skin regularly and effectively – without scheduling appointments or paying extra costs. With 220 light points, the mask covers the entire surface of your face, ensuring effective and even treatment with every session.


Exclusive Series

blue light
470 nm

red light
630 nm

deep red light
660 nm

near infrared light
850 nm

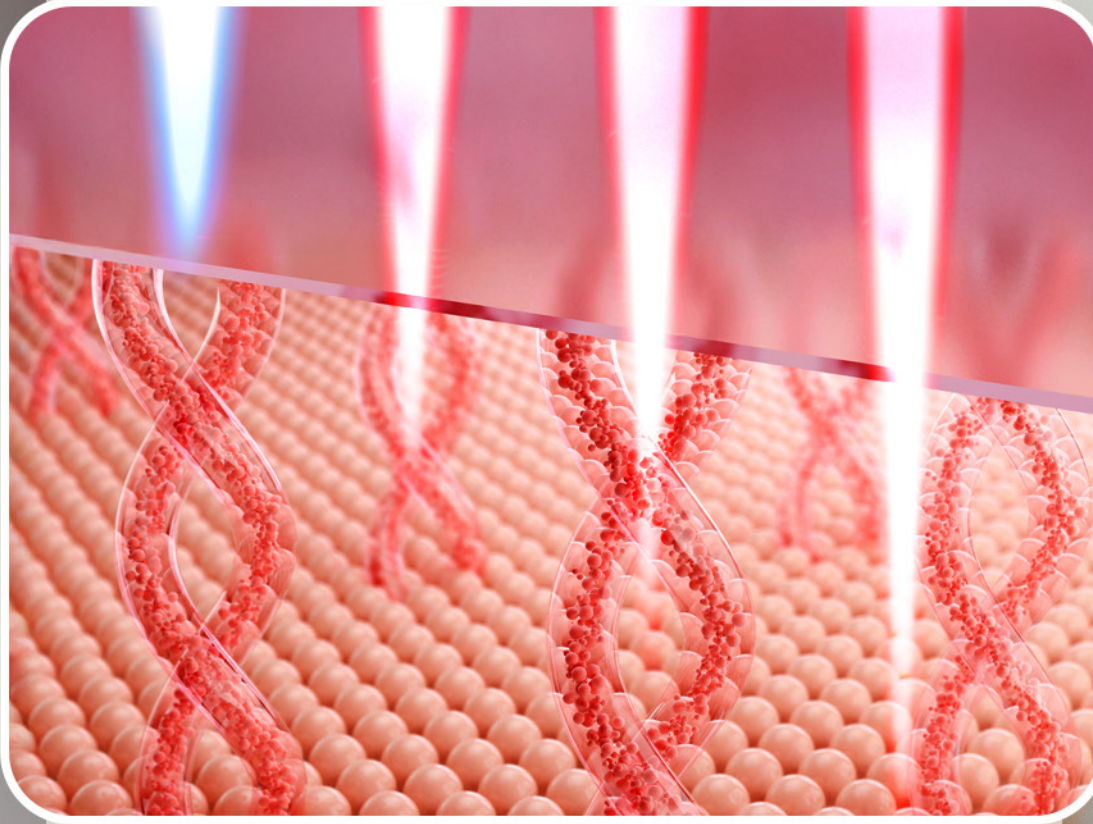

Light adjusted to your skin's needs

Neno Beauty Alva is a modern LED therapy mask that uses four precise wavelengths to comprehensively support your skin. Regular use improves skin's look and condition – increasing elasticity, smoothing fine lines, evening out skin tone, and subtly brightening your complexion. Based on these wavelengths, three special modes have been designed to meet different skincare needs.



Exclusive Series


Skin regeneration boost

SOOTHING MODE

Blue light (470 nm) stimulates skin renewal processes, supporting its natural regeneration. Regular use of this mode improves overall skin condition and has strong antibacterial effects. It reduces inflammation, soothes redness and irritation, and helps regulate sebum production. A great choice when you want to speed up natural repair processes.



★ blue light

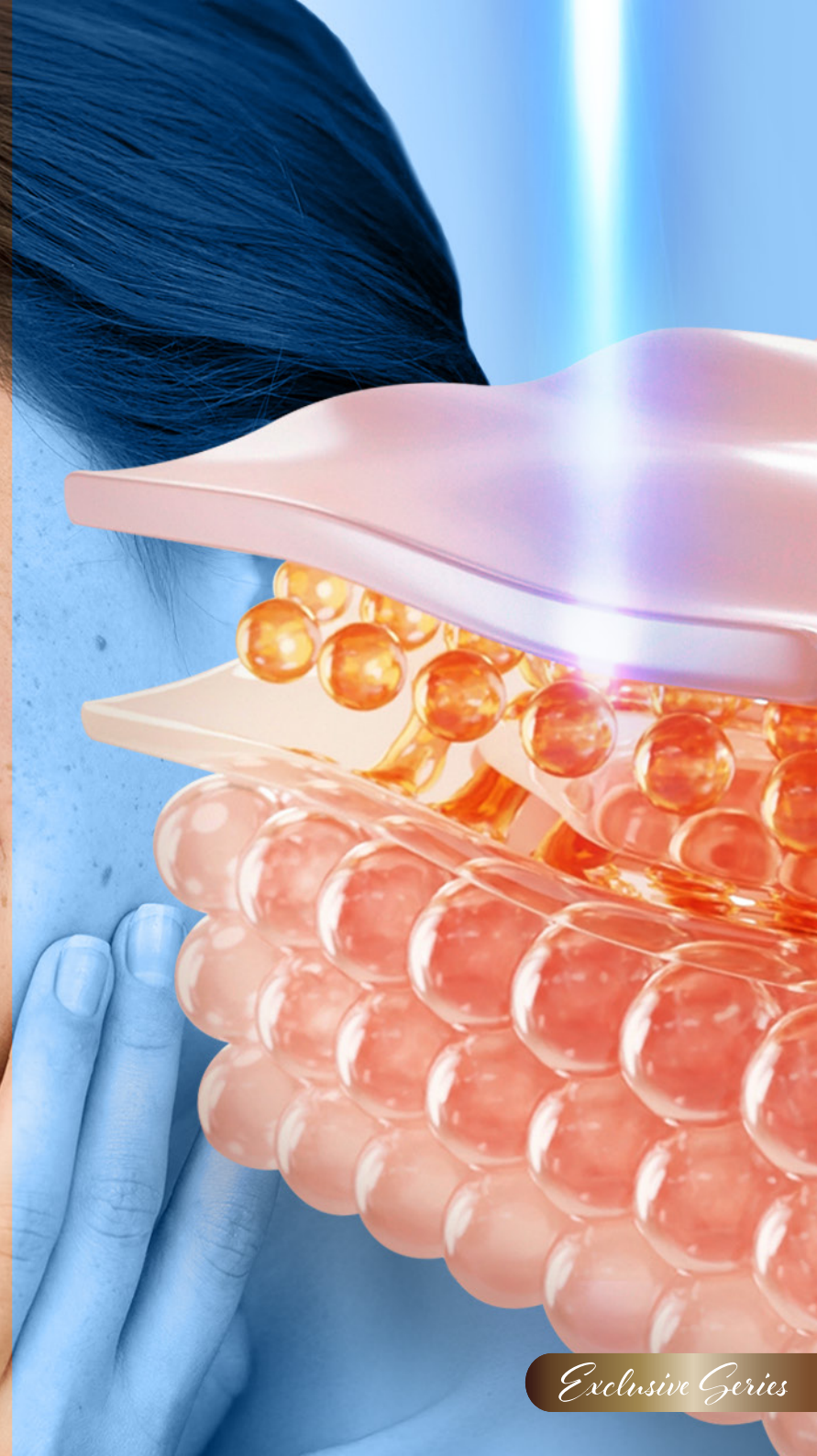

Exclusive Series


Brightening mode

A combination of red light (630 nm) and deep red (660 nm) gives the skin a healthy, radiant glow. This mode stimulates collagen production, supports skin renewal, evens out skin tone, and helps reduce discoloration. Your skin becomes firmer and more luminous.



✦ red light

✦ deep red light

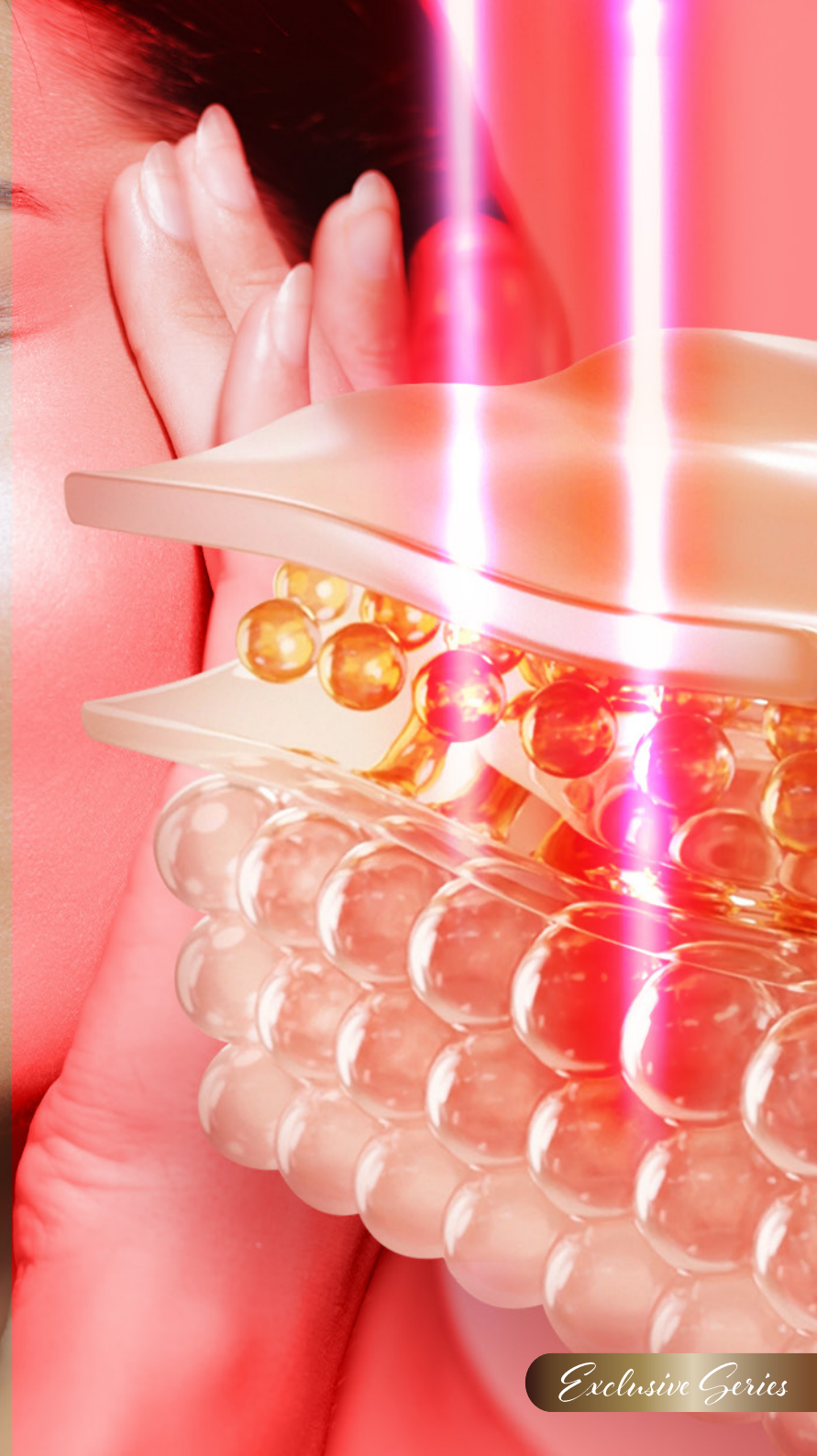

Exclusive Series

Anti-aging mode

The joint action of red light (630 nm), deep red (660 nm), and near-infrared (850 nm) penetrates deeper layers of the skin. This program improves firmness, reduces fine wrinkles, and helps slow down the aging process – leaving skin looking younger and healthier.

✦ red light

✦ deep red light

✦ near infrared light

Exclusive Series



220 LEDs

The Neno Beauty Alva mask is equipped with 220 LED lights. Thanks to this high number of light points, the device covers the entire face surface, ensuring even and intensive action. Each session is fully effective, delivering just the right amount of light your skin needs to regain a healthy, radiant look.

Exclusive Series



neno
beauty

neno
beauty



Timer

The wireless design allows you to use the mask however you like – while reading or relaxing on the couch. The automatic timer ensures the correct treatment time and signals the end of each 10-minute session, so you can fully enjoy your relaxation.

Exclusive Series

Cordless

Thanks to its wireless design and the head strap supporting the mask, you can take care of your skin at any moment – while resting, reading, or doing household chores. Easy charging with a USB-C cable ensures the mask is always ready to use.



Exclusive Series

Intuitive operation

With just one press of a button, the treatment begins. Voice prompts inform you about the selected program, so you always know you're using the right mode. Simple, convenient, and hassle-free.



Exclusive Series

Comfortable design

Neno Beauty Alva is designed for both comfort and safety. Around the eyes, soft silicone gently adheres to the skin, while the flexible strap helps adjust the mask to your head shape. A chin support ensures stability so the session runs smoothly.



High-quality silicone eye protection

Chin support



Flexible, adjustable strap holding the mask on the face

Exclusive Series

Safety guidelines

Prepare your skin:

Before using the LED mask, thoroughly wash and dry your face. Use on clean skin without applying skincare products.

After intensive treatments:

Allow a few days' break between LED mask use and strong treatments like retinoids or chemical peels. The length of the break depends on product strength – when in doubt, consult a cosmetologist.

Photosensitivity warning:

Do not use the LED mask if you are taking photosensitizing medications (e.g. certain antibiotics, retinoids) or herbs that increase light sensitivity (e.g. St. John's wort, chamomile, calendula, angelica).

Pregnancy and breastfeeding:

Use of the LED mask during pregnancy is not recommended. Consult a doctor before use. During breastfeeding, as skin may be more sensitive, a medical consultation is also recommended.

Protect your eyes and follow session time:

Do not look directly at the LEDs, always use the mask with intact eye protection, and never exceed the recommended session time.

Usage frequency:

Use the LED mask 3–5 times per week (for 10–20 minutes), choosing the mode that suits your needs. Before use read all safety precautions included in the user manual.





Hygienic storage

After each session, wipe the inside and outside of the mask with a dry or slightly damp cloth, then store it in the included pouch. This convenient solution protects the device from dust and scratches while keeping it clean for the next use.

neno:
beauty

Exclusive Series

neno[®] beauty

ALVA



SPECIFICATION

Power supply:	5V 0.8A
Power:	6W
Battery:	1000mAh
Timer:	10 min
Material:	ABS+PC+silicone
Weight:	250 g
Dimensions:	187x177x119 mm
Box dimensions:	230x230x123 mm
Box weight:	1.05 kg
EAN code:	5902479677156
Set contents:	LED mask, storage pouch, type-C cable, user manual
Contact:	tel: 731 594 301; mail: info@neno.pl
Producer:	KGK Trend sp. z o.o. Ujastek 5B, 31-752 Cracow

AWARDS FOR
KGK TREND COMPANY
Owner of the Neno brand:



If you don't know which product to choose, need help or support, or have questions for our experts, call or write directly to the Neno Customer Support Center. Helpline: 731 594 301, mail: info@neno.pl