

neno[®] beauty

LUVE

MULTIFUNCTIONAL FACIAL CARE DEVICE
WITH MICROCURRENTS
AND SKIN REJUVENATION TECHNOLOGY

Exclusive Series




Exclusive Series

Discover the Neno Beauty series

Neno Beauty combines modern technology with the comfort of everyday care. The devices in this series make skincare simple, effective, and tailored to your lifestyle. Without visiting a salon, you can enjoy the effects of advanced treatments and rediscover your skin's natural glow.




Exclusive Series

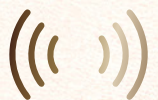

AIR SHOT
mode



Skin contact
sensor

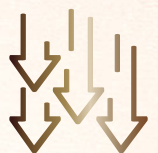


BOOSTER
mode



5 levels
of intensity



MC
mode

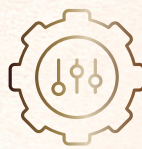

4 modes



DERMA SHOT
mode



Timer



Exclusive Series

Professional Skincare at Home

Neno Beauty Luve is an advanced facial care device that combines 4 modes and 5 intensity levels, giving you full control over every treatment. Discover what Neno Luve can do for you:

AIR SHOT

Blue LED light (460 nm), electrical micro-impulses, and short-wave vibrations improve skin texture, give it a fresh look, and support better absorption of active ingredients in later care steps.

BOOSTER

Yellow LED light (590 nm), gentle vibrations, and electroporation help active ingredients penetrate deeper, enhance skincare effectiveness, improve microcirculation, and leave skin hydrated and radiant.

MC

Green LED light (521 nm), vibrations, and microcurrents stimulate collagen production, improve skin firmness, and smooth fine lines and wrinkles, restoring a youthful look.

DERMA SHOT

Red LED light (630 nm), vibrations, and EMS technology firm facial muscles, enhance skin tone, and refine facial contours.



Exclusive Series

Electrical micro-impulses

Neno Beauty Luve uses electrical micro-impulses that work at the cellular level. Gentle impulses stimulate microcirculation, improve oxygenation, and support regeneration processes. The skin becomes nourished, radiant, and healthy-looking.



Exclusive Series




Electroporation

Electroporation is an advanced, non-invasive technology that uses short electrical pulses to increase the permeability of skin cell membranes. This allows active ingredients in cosmetics to penetrate deeper layers, significantly enhancing skincare effectiveness. With Neno Beauty Luv, electroporation is now available not only in salons but also in your home routine.

Exclusive Series

Microcurrents


Microcurrents in Neno Beauty Luve gently stimulate the skin, supporting microcirculation and oxygenation of tissues. Skin regains firmness, a natural glow, and a healthy look. They also help boost collagen and elastin production, smoothing fine wrinkles and improving tone and condition.

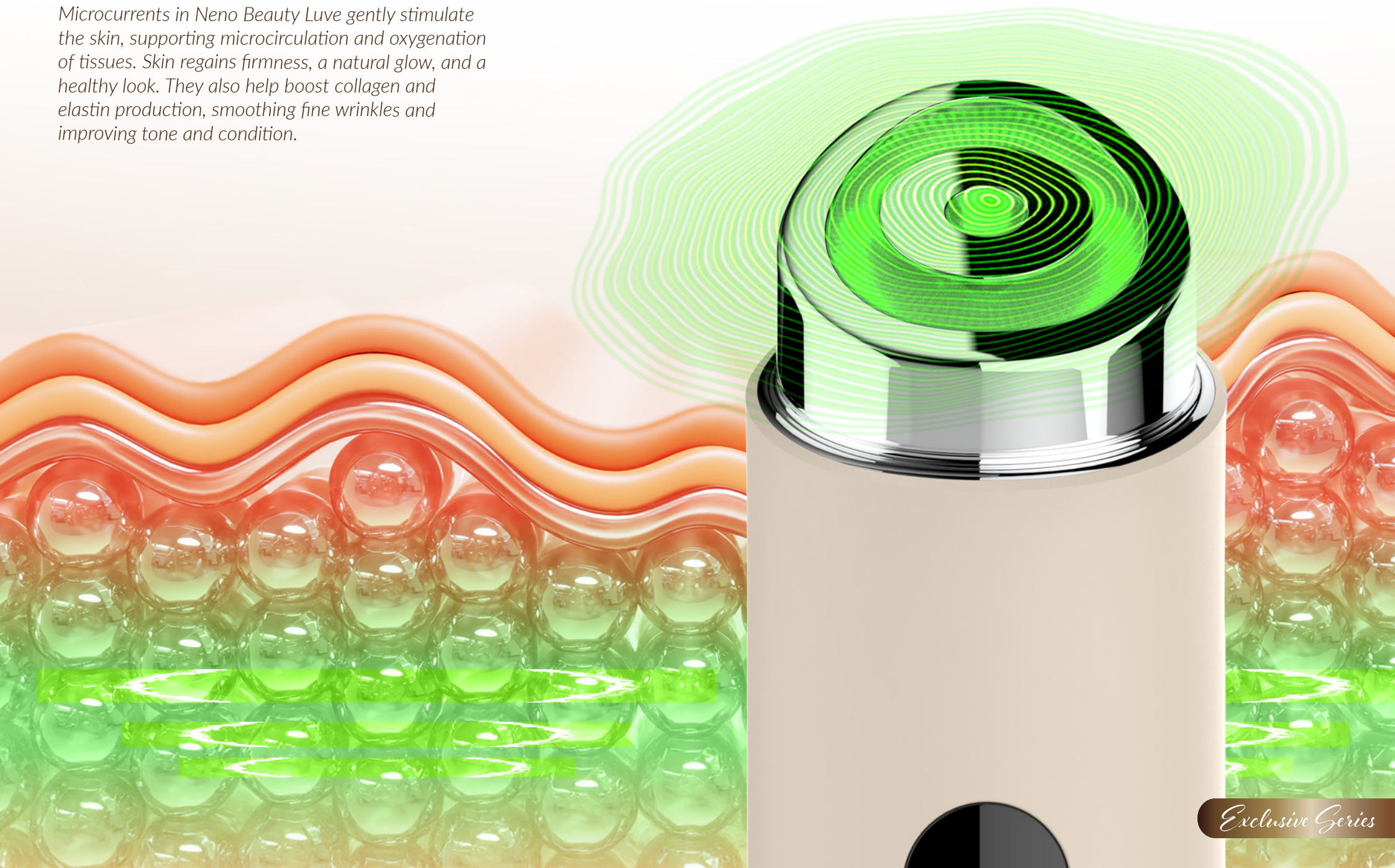

Exclusive Series

EMS technology

EMS is the most intensive form of microcurrent stimulation used in Neno Beauty Luve. It works like a workout for facial muscles, helping maintain their natural tone and supporting skin elasticity. Regular use delivers the best results, depending on skin predispositions. EMS supports firmness, elasticity, and improved facial contours.

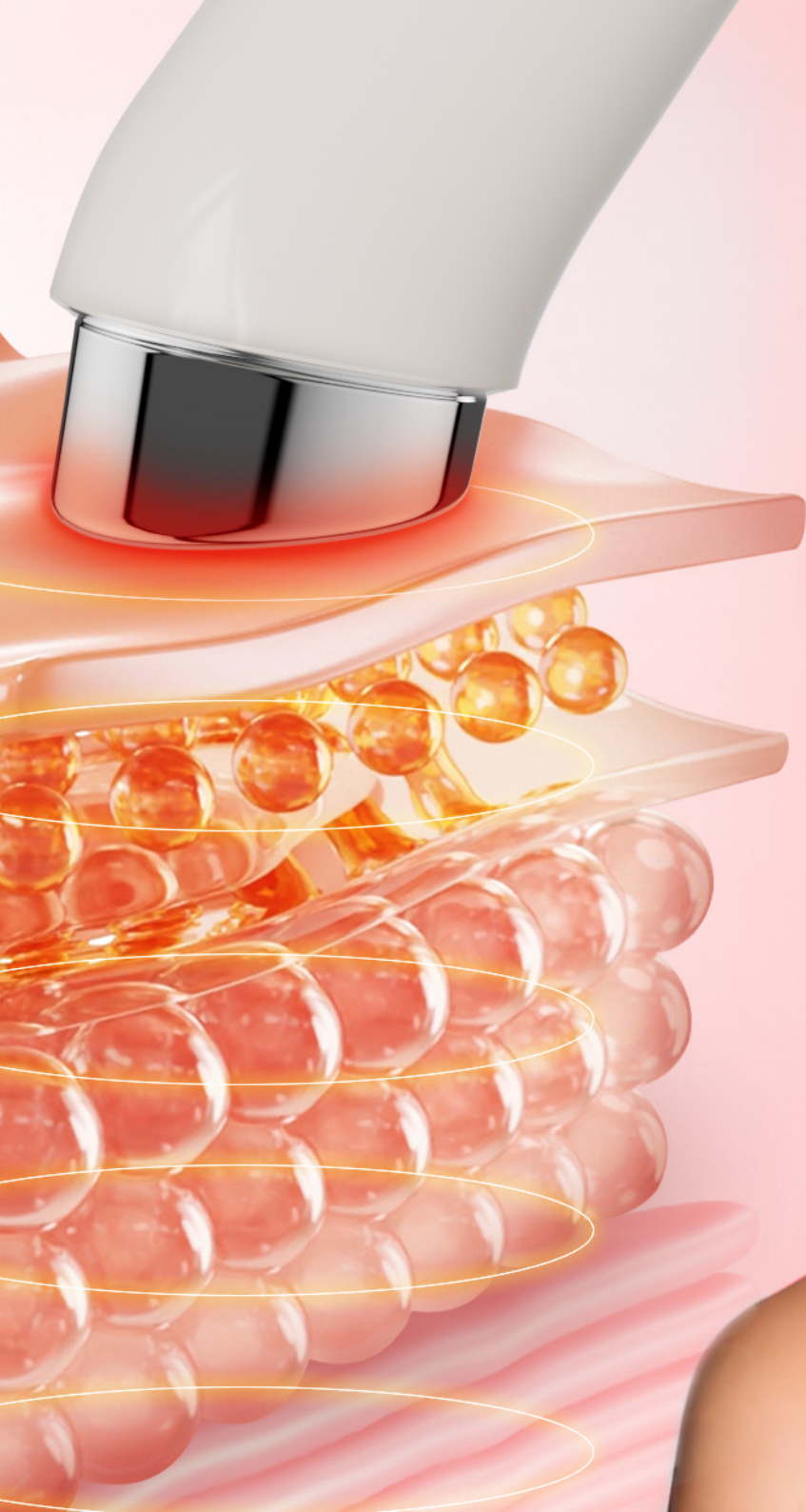

Exclusive Series

Four light wavelengths

460 nm

BLUE LIGHT

supports natural skin repair processes, antibacterial, reduces inflammation, soothes redness and irritation, helps regulate sebum.

521 nm

GREEN LIGHT

reduces redness, soothes the complexion, and improves overall condition.

590 nm

YELLOW LIGHT

calming effect, evens skin tone, gives skin a fresh glow.

630 nm

RED LIGHT

stimulates microcirculation, evens skin tone, reduces discoloration, soothes the complexion.


Exclusive Series

Sonic vibrations

Neno Beauty Luve uses sonic vibrations tailored to each skincare mode, ensuring effective and comprehensive support for daily care. Special vibration patterns enhance ingredient absorption, improve facial contouring, and boost firmness and elasticity.



Exclusive Series

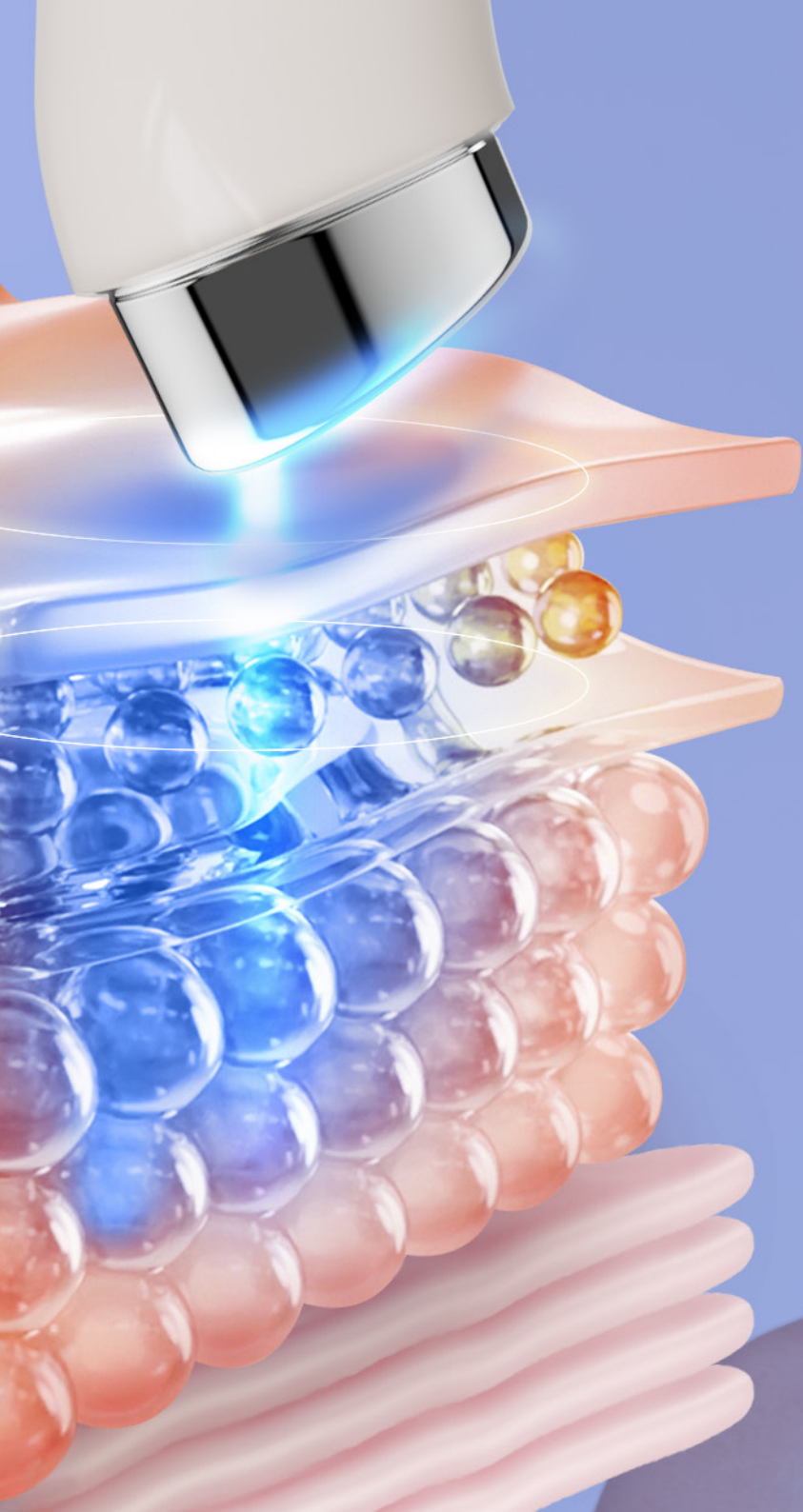

Air Shot

Combines blue LED light, micro-impulses, and short-wave vibrations. Reduces imperfections and improves skin texture. May cause slight tingling or a faint smell of burning tiny hairs – normal and safe.



Exclusive Series



Booster

Combines yellow LED light, vibrations, and electroporation to support deeper absorption of active ingredients. Enhances skincare effectiveness, improves circulation, and evens skin tone. Use with gel-based or water-based cosmetics rich in active ingredients (hyaluronic acid, peptides, niacinamide, ceramides). Tingling may occur – a natural effect.

Exclusive Series

MC

Combines green LED light, vibrations, and microcurrents. Stimulates collagen, improves firmness, and reduces fine lines. Tingling may occur – a normal response to stimulation.



Exclusive Series

Derma Shot

Combines red LED light, vibrations, and EMS. Firms muscles, tones skin, and enhances facial contours.

May cause gentle twitches or tingling – a natural result of impulses.



Exclusive Series



Five intensity levels

Each mode offers 5 intensity levels, letting you personalize treatments. Start at the lowest level and gradually increase over weeks of use. For delicate areas (eyes, nose, lips), always use the lowest level.

Exclusive Series

Massage directions

AIR SHOT – Use the edge or tip of the head, not the whole surface. Move slowly along chosen areas or apply to specific spots. Use on clean, dry skin without products.

BOOSTER – Apply the full head surface in circular motions (from the center of the face outward and upward). Always use with light gel-based skincare rich in active ingredients. Add more product if it absorbs fully during use.

MC – Apply the full head surface, moving slowly upward and outward. Use with oil-free conductive gel for EMS or light, oil-free skincare. Reapply gel if it dries.

DERMA SHOT – Apply the full head surface, moving slowly upward and outward. Pause for 1–2 seconds on selected spots. Always use with oil-free conductive gel or light, oil-free skincare. For best results, use MC mode before Derma Shot.



Treatment order

The device allows use of all four modes in one cycle.
For best results and safety, follow this order: Air Shot
→ Booster → MC → Derma Shot. Only Air Shot is
always used on clean, dry skin without products.



Exclusive Series

Safety and hygiene

- Clean after every use – wipe with a dry or slightly damp cloth with mild soap. For disinfection, use isopropyl alcohol wipes (up to 70%). Avoid harsh detergents, bleach, or immersion in water.
- Do not use if you have electronic implants (e.g. pacemaker, insulin pump) or metal implants in the head/neck (dental, braces, hearing implants).
- Avoid sensitive areas – do not use directly under eyes, on eyelids, or thyroid area. Around eyes, nose, mouth – always use the lowest intensity.
- Consistency – best results with regular use:
Air Shot: 2–3 times per week (max 5 min)
Booster: 1–3 times daily (5–10 min)
MC: 1–3 times daily (5–10 min)
Derma Shot: 1–3 times daily (5–10 min)
- Do not use during pregnancy or breastfeeding.
- Always read all safety precautions in the manual before use.



Exclusive Series

Skin contact sensor + Timer

Neno Beauty Luve has a built-in 5-minute timer and a skin contact sensor. The timer starts only when the head touches the skin, and LED lights switch on only on contact – ensuring safety and avoiding accidental light exposure to the eyes.

5 min



Exclusive Series

Accessories included

Neno Beauty Luve comes with a practical protective cap for the massage head and an elegant pouch for safe storage. The cap keeps the device clean and ready to use, while the pouch is perfect at home or on the go.



Exclusive Series

neno[®] beauty LUVÉ

TECHNICAL DATA

Model:	Neno Beauty Luve
Power supply:	5V 0.7A
Power:	5W
Battery:	2200 mAh
Timer:	5 min
Material:	ABS+PP+304 stainless steel
Dimensions:	38.2×182.7×38.6 mm
Weight:	304g
EAN code:	5902479677170
Box dimensions:	148×245×60 mm
Box weight:	0,64 kg
PN:	NEN-DOM-PI003
Set contents:	massager, storage pouch, microfibre cloth type-C cable, user manual
Contact:	Tel.: 731 594 301; mail: info@neno.pl
Producer:	KGK Trend sp. z o.o. ul. Ujastek 5B, 31-752 Cracow



AWARDS FOR
KGK TREND COMPANY
Owner of the Neno brand:



If you do not know which product to choose, need help or support, or have questions for our experts, please call or write directly to the Neno Customer Support Centre.
Hotline: 731 594 301, mail: info@neno.pl