nenò[®]

Pulito

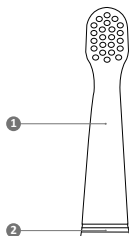
Table of contents


Spis treści

4	Instrukcja obsługi	PL
6	User manual	EN
8	Bedienungsanleitung	DE
11	Návod k použití	CZ
13	Használati utasítás	HU
16	Používateľská príručka	SK
18	Användarhandbok	SE
20	Käyttäjän käsikirja	FI
23	Brukerhåndbok	NO
25	Brugermanual	DK
27	Gebruikershandleiding	NL
30	Manual del usuario	ES
32	Manuale utente	IT
35	Manual de utilizare	RO
37	Korisnički priručnik	HR/BA
39	Uputstvo za upotrebu	RS/ME
42	Uporabniški pravilnik	SI
44	Εγχειρίδιο χρήστη	GR
47	Lietotāja rokasgrāmata	LV
49	Vartotojo vadovas	LT
51	Kasutusjuhend	ET
54	Посібник користувача	UA
56	Ръководство за потребителя	BG
59	Упатство за употреба	MK
61	Manual do utilizador	PT
63	Manuel de l'utilisateur	FR


1


2

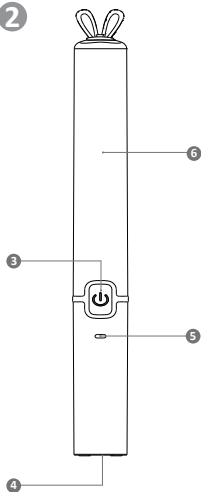3

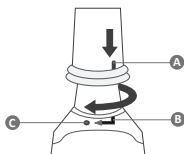4

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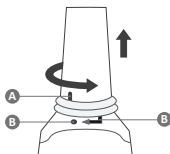EN

USER MANUAL

Dear customer, thank you for purchasing NENO PULITO. The product you have purchased is an electric toothbrush for children. This device is used for brushing teeth and maintaining oral hygiene. Please read the instructions below before using the device.

01. PRECAUTIONS

1. If you do not understand the instructions for use or the precautions listed in this manual, do not use the appliance. Incorrect use of the appliance may damage it.
2. Do not take the product apart. This may result in damage to the product, electric shock and invalidation of the warranty.
3. Take care that the device does not fall to the ground. If it has been damaged, contact the manufacturer.
4. Do not use accessories other than those recommended by the manufacturer.
5. Do not use flammable or explosive substances to clean the product.
6. This product is for brushing teeth only. Do not use the appliance for any other purpose.
7. Discontinue use of the appliance if it is visibly damaged.
8. Do not place the product in water. This may lead to a malfunction and risks fire, electric shock or explosion.
9. Keep the appliance away from containers with water, such as a bathtub or washbasin.
10. If the toothpaste you use contains peroxides or bicarbonates (often present in tooth whitening toothpaste), wash the appliance thoroughly after use to avoid damaging its plastic casing.
11. Protect the device's charging port from moisture ingress.

12. If the use of the device causes discomfort to your child or bleeding gums, consult your dentist before further use.
13. The appliance is not intended for use by persons with limited physical, cognitive or mental capabilities or who lack experience or knowledge, including children. If such a person wishes to use the appliance, they must be supervised or instructed in the instructions for use by a person responsible for them.
14. The appliance is not a toy.
15. Do not expose the toothbrush to direct sunlight.
16. The device is for personal hygiene and should not be used by more than one person.

02. DESCRIPTION OF THE DEVICE

FIG. 1

1. Head with toothbrush
2. Rubber gasket

FIG. 2

3. Powering the toothbrush
4. Charging the toothbrush
5. LED
6. Closing the toothbrush

03. CONTENTS OF THE KIT

1. Neno Pulito toothbrush
2. Extra brush head
3. Charging cable
4. Instructions for use

04. HOW TO USE

1. Pull off the toothbrush cap (FIG.3-6)
2. Place the toothbrush head on the holder, if not already fitted.

NOTE:

- To fit the brush head, slide element A on the head onto element B of the brush holder and twist it sideways until it sits on element C (FIG. 3).
 - To remove the toothbrush head, twist element A on the head sideways until it is on element B and pull it upwards (FIG. 4).
3. Squeeze the toothpaste onto the toothbrush.
 4. Press the power button (FIG.2 - 3). By pressing the button, select one of the following modes (LED FIG.2 - 5, will change colour):
 - Green colour - gentle mode - used for cleaning sensitive teeth and tongue.
 - Red colour - standard mode - used for cleaning teeth.
 - Blue colour - massage mode - for gum stimulation.
 5. Gently press the toothbrush against your teeth and gums..

NOTE: Thoroughly clean every part of your child's teeth. Every 30 seconds of brushing, the toothbrush stops for 1 second to allow the user to change the area being brushed. Correct tooth brushing should take a total of 2 minutes. The device will automatically switch off after 2 minutes.

6. When you have finished brushing your teeth, switch off the toothbrush by pressing the power button again (FIG.2 - 3).
7. After brushing your teeth, wash the toothbrush under running water. Shake the remaining water out of the brush bristles and use a soft cloth to wipe the handle of the unit.

05. CHARGING YOUR TOOTHBRUSH

1. Plug the cable from the kit into the port located at the base of the toothbrush (FIG.2 - 4).
2. The lamp will illuminate red while charging.
3. When the unit charges, the light will turn green.

06. CLEANING THE TOOTHBRUSH

1. Wash the brush head and bristles after using the device. Allow them to dry naturally. Do not close the toothbrush.
2. Do not push or pull on the bristles with your fingers. This can deform the toothbrush and break it.
3. Wipe the handle with a soft material and let it dry naturally.
4. Do not submerge the handle under water.
5. Do not use water above 40°C or acidic cleaners to clean the device.

07. SPECIFICATIONS

Power supply: DC1.2V

Battery: 700mAh

Waterproof: IPX7

Brush bristles: Nylon

Vibration per minute: 38000

Built-in timer: YES

Dedicated age: 3+

Dimensions: 20x19x193mm

Weight: 73g

Charging time: up to 6 hours

08. WARRANTY CARD

The product comes with a 24-month warranty. Warranty terms and conditions can be found at: <https://neno.pl/gwarancja>

Details, contact and service address can be found at: <https://neno.pl/kontakt>

Specifications and contents are subject to change without notice. We apologise for any inconvenience.

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EN

The crossed out trash can symbol indicates that unusable electrical or electronic devices, its accessories (such as power supplies, cords) or components (for example batteries, if included) cannot be disposed of alongside with household waste. In order to dispose of the devices or its components (for example, batteries) deliver the device to the collection point, where it will be accepted free of charge. Disposal is subject to the recast version of the WEEE Directive (2012/19/ EU) and the Directive on batteries and accumulators (2006/66 / EC). Proper disposal of the device prevents degradation of the natural environment. Information about the collection points of the facilities is issued by the competent local authorities. Incorrect disposal of waste is subject to penalties provided for by the law in force in the given area.

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