**nenor**

**Pulsar**


## Table of contents


Spis treści

|    |                            |       |
|----|----------------------------|-------|
| 4  | Instrukcja Obsługi         | PL    |
| 5  | User Manual                | EN    |
| 7  | Bedienungsanleitung        | DE    |
| 8  | Návod K Použití            | CZ    |
| 9  | Használati Utasítás        | HU    |
| 11 | Používateľská Príručka     | SK    |
| 12 | Användarhandbok            | SE    |
| 14 | Käyttäjän Käsikirja        | FI    |
| 15 | Brukerhåndbok              | NO    |
| 17 | Brugermanual               | DK    |
| 18 | Gebruikershandleiding      | NL    |
| 19 | Manual Del Usuario         | ES    |
| 21 | Manuale Utente             | IT    |
| 22 | Manuel De L'utilisateur    | FR    |
| 24 | Manual De Utilizare        | RO    |
| 25 | Korisnički Priručnik       | HR/BA |
| 27 | Uputstvo Za Upotrebu       | RS/ME |
| 28 | Uporabniški Pravilnik      | SI    |
| 30 | Εγχειρίδιο Χρήστη          | GR    |
| 31 | Lietotāja Rokasgrāmata     | LV    |
| 32 | Vartotojo Vadovas          | LT    |
| 34 | Kasutusjuhend              | ET    |
| 35 | Посібник Користувача       | UA    |
| 37 | Ръководство За Потребителя | BG    |
| 38 | Упатство За Употреба       | MK    |
| 40 | Manual Do Utilizador       | PT    |


1

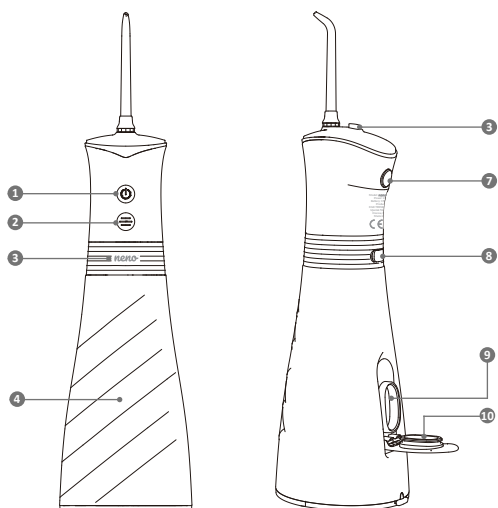
2

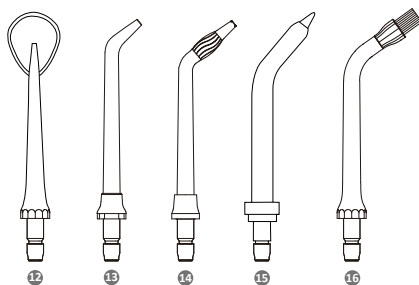3

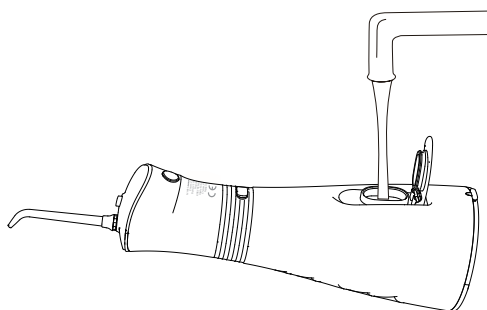4

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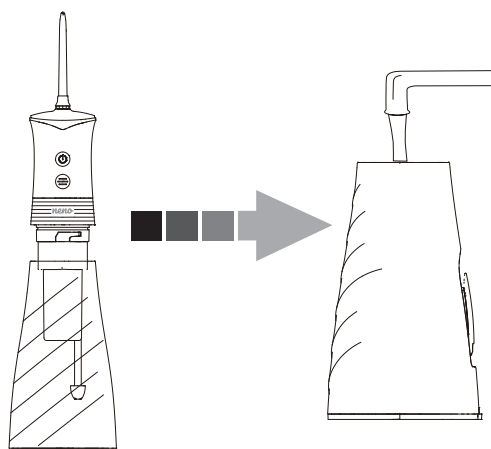# EN

## USER MANUAL

Dear Customer, thank you for purchasing the Neno Pulsar. This product is an electric oral irrigator designed for adults. This device is used to maintain oral hygiene. Please read these instructions carefully before use.

## 01. PRECAUTIONS

1. Do not touch the charging port or power adapter plug with wet hands. Do not connect or disconnect the power adapter from the electrical outlet with wet hands.
2. Ensure that the charging port is dry before charging. Always use the original charging cable.
3. The product is not intended for use by persons with physical, sensory or mental limitations or lack of experience, such as children, unless they are under appropriate supervision. Caregivers should ensure the user's safety and that children do not play with the device.
4. The product contains a battery. Do not throw it into a fire, expose it to heat, or store it at high temperatures or near heat sources.
5. Do not immerse the product in water. Prolonged contact with water may damage the device.
6. Do not disassemble or modify the device, cable or charger. Do not place heavy objects on the cable or pull on it.
7. If the device malfunctions or fails to operate properly, do not use it.
8. If there are any malfunctions or failures, do not use the product. Contact an authorised service centre.
9. Always disconnect the charger and cable from the power supply before cleaning the device.
10. Do not use the device without the tip installed.
11. Keep the device out of the reach of children and infants. Prevent them from coming into contact with the product.
12. People with oral or gum diseases should consult a doctor before using the device.
13. The product is intended for oral hygiene only. Do not use it for other purposes.
14. Do not use the product at temperatures above 40°C to avoid damage.
15. Do not share the tips with other people.
16. Only use original tips recommended by the manufacturer.
17. When disconnecting the power supply, hold the plug – do not pull on the cable.
18. Clean the device regularly according to the instructions.
19. If the tip becomes worn or damaged, replace it immediately.
20. Do not open, disassemble or attempt to repair the device yourself.
21. Do not fill the tank with sea water, brine or concentrated essential oils that do not dissolve in water. Using such liquids may damage the device. If you use mouthwash, dilute it appropriately before use.

## 02. DESCRIPTION OF THE DEVICE

### Figure 1

1. Device on/off button.
2. Mode selection button.
3. Mode/charging indicator.
4. Water tank.
5. Irrigator nozzle.
6. Nozzle release button and socket.
7. Manual mode/pause button.
8. Charging port.
9. Water filling/emptying opening.
10. Opening closure.

Figure 2

- 11. Tongue cleaning nozzle – suitable for removing tongue coating.
- 12. Standard nozzle – suitable for daily cleaning.
- 13. Tornado nozzle – rotating stream for better cleaning.
- 14. Periodontal pocket cleaning nozzle - for people with periodontitis or other gum problems.
- 15. Orthodontic nozzle – for people with braces.

03. SET CONTENTS

- Neno Pulsar irrigator
- Replaceable nozzles – 5 pieces
- Fabric travel bag
- USB-C cable
- User manual

04. HOW TO USE

1. Insert one of the nozzles included in the set into the socket located at the top of the irrigator by pressing the nozzle release button (Figures 1-6).
2. Fill the tank with water using one of the following two methods:
  - Open the cap (Figure 1-10) on the water tank (Figure 1-4) and, holding the device horizontally, pour water through the opening (Figure 1-9) as shown in Figure 3.
  - Remove the water tank (Figure 1-4) by gently turning it to the side. Holding it vertically, fill it with water as shown in Figure 4.
3. Start the device using the power button (Figure 1-1). The mode indicator (Figure 1-3) will light up and you will be able to select one of the available modes using the selection button (Figure 1-3):
  - Clean - continuous water flow - continuous white light. Duration: two minutes.
  - Soft - continuous water flow (gentler) - continuous blue light. Duration: two minutes.
  - Massage – 3 water jets per second (gum massage) – flashing blue light. Duration: two minutes.
  - DIY - continuous stream with adjustable power - steady light blue light.

**NOTE:** The device has a mode memory function.

**NOTE:** During operation, you can pause the device by pressing and holding the pause button. (Figure 1-7). When the button is released, the device will resume operation. The pause time is not included in the operation mode timer.

**NOTE:** In DIY mode, you can manually adjust the spray power by holding down the mode button (Figure 1-3). When you release the button, the setting will be automatically saved.

**NOTE:** The button (Figure 1-7) is used for manual operation of the device. Pressing the button starts continuous spraying, and releasing it stops spraying, allowing you to adjust the duration of the stream.

4. After use, press the off button (Figure 1-1) to turn off the device.

**NOTE:** The device will switch off automatically after 2 minutes.

5. Open the tank opening cap (Figure 1-10) and pour out any remaining water.
6. Press the on button (Figure 1-1) to pump the water out of the pipes.
7. Switch off the appliance again, close the tank opening (Figure 1-10) and wipe the appliance dry.

05. CHARGING

1. Insert the charging cable into the charging port on the device (Figure 1-8).
2. The charging indicator (Figure 1-3) will display one of the following colours:
  - **Red LED (steady light)** - charging in progress.
  - **Green LED** - charging complete.
  - **Red LED (flashing)**: low battery level.
  - **Red LED (fast flashing)**: battery overheating – stop charging immediately.

**Note:** A full charge takes approximately 6 hours. Once fully charged, the device automatically stops charging to prevent overcharging.

06. CLEANING

Cleaning individual components:

**Device body**

- Use a dry or slightly damp cloth.
- Do not immerse the body in water.

**Water tank**

- Wipe with a dry or damp cloth, or remove and rinse with clean water.
- Remove water from the tank.

**Tips**

- Clean with a dry or damp cloth, or remove the nozzle and rinse it with water, then dry.
- Do not bend, pull or twist the nozzle.
- Recommendation: replace the tip every 6 months.

07. TROUBLESHOOTING

| PROBLEM                             | POSSIBLE CAUSE                                       | SOLUTION                                                |
|-------------------------------------|------------------------------------------------------|---------------------------------------------------------|
| The device does not work            | Newly purchased or unused for more than 3 months     | Charge the battery to 100%                              |
| Low water pressure                  | Blocked or deformed tip                              | Blocked or deformed tip                                 |
| Short operating time after charging | Worn-out battery<br>Charging time too short          | Contact customer service<br>Charge for at least 4 hours |
| Charging indicator not working      | Incorrect power connection<br>Charging circuit fault | Check the power connection<br>Contact customer service  |
| No water flow                       | The appliance has not been used for a long time      | Remove the nozzle and rinse the nozzle hole             |

08. SPECIFICATIONS

**Power supply:** 5V - 0,5A

**Power:** 2.5W

**Dimensions:** 80x86x230mm

**Weight:** 308 g

**Charging time:** approx. 6 hours

**Battery:** 1500mAh, lithium-ion

**Material:** ABS, PMMA

**Charging port:** USB-C

**Water resistance:** IPX7

09. WARRANTY CARD

This product is covered by a 24-month warranty. The warranty terms and conditions are available at:

<https://nenopl/gwarancja>

For details, contact information, and service addresses are available at: <https://nenopl/kontakt>

Specifications and kit contents are subject to change without prior notice. We apologise for any inconvenience.

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