nenor[®]

Integro

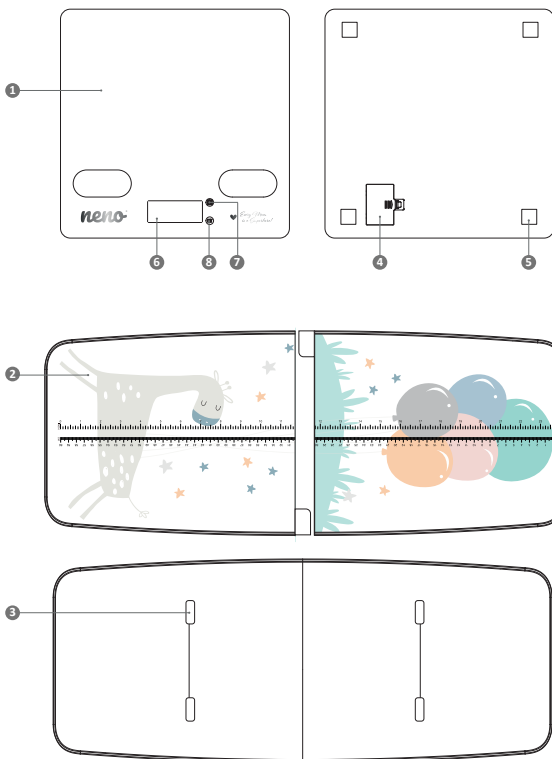
Table of contents


Spis treści

5	Instrukcja obsługi	PL
9	User manual	EN
12	Bedienungsanleitung	DE
16	Manuel de l'utilisateur	FR
20	Návod k použití	CZ
23	Používateľská príručka	SK
27	Használati utasítás	HU
30	Användarhandbok	SE
34	Käyttäjän käsikirja	FI
37	Brugerhåndbok	NO
41	Brugermanual	DK
44	Gebruikershandleiding	NL
48	Manual del usuario	ES
51	Manuale utente	IT
55	Manual de utilizare	RO
59	Korisnički priručnik	HR/BA
62	Uputstvo za upotrebu	RS/ME
66	Uporabniški priročnik	SI
69	Εγχειρίδιο χρήστη	GR
73	Lietotāja rokasgrāmata	LV
77	Vartotojo vadovas	LT
80	Kasutusjuhend	ET
84	Керівництво користувача	UA
87	Ръководство за употреба	BG
91	Упатство за употреба	MK
95	Manual do utilizador	PT

A


B

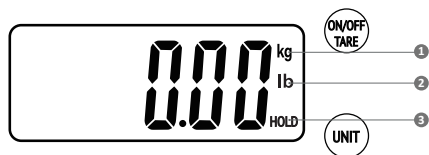
C

1

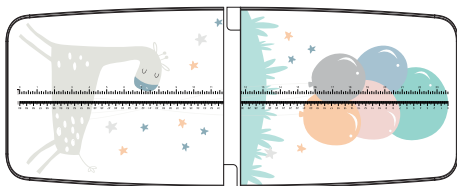2

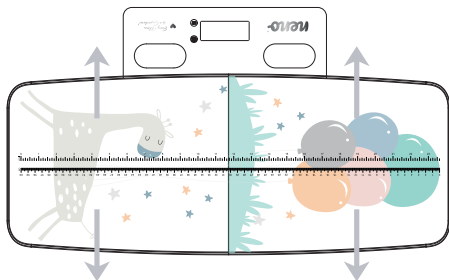D

1

2

3

4

3 s.

5

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Szanowny Kliencie,

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USER MANUAL

Dear Customer,

Thank you for purchasing the Neno Integro 3-in-1 electronic bathroom scale with app and baby weighing pad. The product is designed to allow accurate and convenient weight measurement of infants, children and adults. In combination with a dedicated mobile application, the device also allows you to monitor body parameters such as BMI, percentage of body fat, muscle tissue, water in the body and others. This user manual contains all the information necessary for proper use. Please read it and keep it if you need to use it again.

01. SAFETY PRECAUTIONS

1. Use the device only for its intended purpose, i.e. for measuring body weight.
2. Before first use, read the instruction manual carefully and keep it for later reference.
3. The product is not a toy – children can only use it under adult supervision.
4. Use only the original accessories that come with the device.
5. Store the scale in a dry place, away from moisture and out of the reach of children.
6. Avoid storing the appliance in places exposed to extreme temperatures (e.g. direct sunlight, near radiators).
7. Do not use the scale in an environment with increased humidity or immediately after bathing – this may affect the operation of electronic components.
8. Do not immerse the device in water or other liquids.
9. Only set the scale on a stable, hard, and level surface.
10. Do not use the device on soft ground (e.g. carpet) as this may affect the accuracy of the measurement.
11. Make sure the surface of the scale is clean and dry before stepping on it.
12. Do not step on the scale with wet feet or when its surface is damp – this may cause slipping.
13. Do not stand on the edges of the platform to avoid tilting it.
14. Elderly people, with limited mobility or disabilities should use the device with the support of another person and, if possible, with additional support.
15. Do not exceed the maximum allowable load specified in the appliance specification.
16. Do not place any objects on the scale when it is not in use.
17. Check the technical condition of the device regularly. If problems occur, stop using the device and contact an authorized service center.
18. Do not drop the product, disassemble the device, do not make repairs or modifications yourself.
19. Use only batteries that are in accordance with the manufacturer's specifications.
20. Make sure the batteries are installed correctly.
21. Check the condition of the batteries regularly and replace them as soon as they are used up to prevent possible leakage.
22. Do not insert new or partially used batteries into the appliance at the same time. This may damage the device.
23. It is not acceptable to use batteries that are damaged or with visible signs of leakage.
24. If the device will not be used for a long time, remove the batteries from the device to avoid the risk of damage.
25. Do not dispose of either the battery or the product in the mixed municipal waste container. Comply with applicable regulations regarding the disposal of electronic equipment and batteries.
26. Clean the scale and tray with a soft, slightly damp cloth using mild cleaners.
27. After cleaning, make sure the appliance is completely dry before using it again.
28. People with medical implants (e.g. pacemaker) should consult a doctor before using the body composition analysis (bioimpedance) function.
29. The scale with baby attachment is designed for weighing children only.
30. Never leave an infant unattended while measuring.
31. When measuring, make sure that the child is correctly positioned on the surface of the scale.

02. SET CONTENTS

1. Bathroom Scale
2. Baby Cover – 2 Parts
3. Screws x 2
4. Screwdriver
5. User Manual

03. DESCRIPTION OF THE DEVICE (FIG. A)

1. Scale Platform
2. Two parts of the baby pad

3. Splints for inserting the overlay
4. Battery cover
5. Non-slip feet
6. LCD screen
7. Power/Tare Button
8. Unit change/rotate screen button UNIT

DISPLAY (FIG. B)

1. Weight Indication
2. Unit symbol on display (kg/lbs)
3. HOLD symbol – means pausing the measurement

04. INSTALLATION

Before first use, carefully remove all components from the packaging, taking special care not to damage the device or its accessories. After unpacking, check the completeness of the kit and make sure that the device does not have visible mechanical damage, such as cracks, deformation or loose elements. If any abnormalities are found, do not use the appliance and contact the dealer or manufacturer.

Battery Installation

To install the batteries, turn the unit upside down and remove the battery compartment cover (FIG. A.4). Then place the batteries in the compartment, paying attention to the correct alignment of the poles. The device is powered by two CR2032 batteries (batteries are not included).

Once the battery is properly placed, close the battery compartment cover, making sure it is properly secured.

Folding Baby Weighing Tray

The tray consists of two pieces that must be connected to each other before the first use.

First, align both parts of the overlay so that the tab in one part is placed in the corresponding hole of the other part. Once the components are properly fitted, screw them together using the included screws.

Do this symmetrically on both sides of the tray, making sure that all connections are stable and properly tightened (FIG. C.1).

Tray Assembly and Disassembly

To mount the tray on the scale, slide the scale into the guides on either side of the bottom of the tray, moving it along the guides until it is properly seated (FIG. C.2).

To disassemble the tray, slide it out in the opposite direction, gently guiding it along the same guides.

05. DEVICE OPERATION

Inclusion

The device can be switched on by overpressure of the power button FIG. D.1

The scale also starts automatically when a load exceeding 5 kg is detected.

Adult Body Weight Measurement

To ensure correct and accurate measurement, place the balance on a hard, stable, and flat surface.

Then gently stand in the middle of the platform. The device will turn on automatically.

During the measurement, stand still until the result stabilizes. When the measurement is complete, the display indicator flashes three times and then the HOLD function is activated to lock the result on the screen.

Do not leave the balance before you have finished taking a measurement, as this may affect the accuracy of the reading.

Measuring your baby's body weight

Install the infant weighing tray. To do this, slide the balance body into the guides located on the bottom of the tray. Gently slide the device along the guides until it is fully and firmly seated. Make sure the scale has been properly attached and does not move around during use.

After installing the tray, place the whole thing on a stable, hard and flat surface, and then turn on the device with the power button.

Carefully place the baby in the middle of the tray, taking extra care and constant supervision.

When the measurement is complete, the result on the display will flash three times and then the HOLD function will be activated, which locks the indication, preventing it from changing under the influence of the child's movements.

Do not remove your child from the scale before the measurement is complete, as this may affect the accuracy of the reading.

Tare function

If additional items (e.g. blanket, towel or diaper need to be used), place them on the tray first and then press the tare button (FIG. D.2). The display will show „0.00“ to indicate that the weight of the item has been zeroed and will not be included in the measurement.

It is recommended to perform tare before each measurement with the use of additional accessories. This function should also be used when there is another object on the balance that cannot be removed and that could affect the measurement result.

Changing Weight Units

The device allows you to measure in two units: kilograms (kg) and pounds (lb).

To change the unit of measurement, press the UNIT button (FIG. D.3).

Change the direction in which your data is displayed

To rotate the display by 180°, press and hold the UNIT button for approximately 3 seconds (FIG. D.4).

After doing this, the orientation of the data on the screen will be changed.

Turning off the device

The device turns off automatically after approximately 1 minute of inactivity. Until switched off, the result of the last measurement remains visible on the display.

It is also possible to manually turn off the appliance by pressing and holding the power button for approximately 3 seconds (FIG. D.5).

Overload indicator

If the maximum permissible load is exceeded, the message „OL“ will appear on the display.

In this situation, remove the scale immediately to prevent damage.

Low battery indicator

When the battery level is low, the display will show „LO“.

In this case, stop using the device and replace the batteries with new ones that comply with the manufacturer's specifications.

06. MOBILE APPLICATION

The scale works with a dedicated **Neno mobile application**, which allows you to extend the functionality of the device to include analysis and archiving of measurement results. The application allows you to conveniently monitor body parameters and track changes over time, presenting data in a clear form of graphs and summaries.

Installing the app

To start using the app's features, you need to download it to your mobile device from the app store and then install it. To connect the scale to the app, turn on the device with the power button and then launch the Neno app. In the application menu, select the option to add a new device. The app will automatically start searching for available BT devices. From the list that appears, select the NENO INTEGRO model. In case of connection problems, make sure the scale is within BT range.

App features – adult user

For an adult user, the app allows you to measure, analyze and record a number of body parameters. In addition to body weight, it is possible to monitor BMI, body water content, body fat levels, bone mass, and other body composition parameters available in the system. All results are automatically saved in the user's history, allowing for long-term observation of changes and progress control.

Child sub-account

The app also allows you to create a separate profile (sub-account) for your child. This feature is designed to monitor the development of the youngest users. In child mode, only body weight is recorded. Individual notes can be added to each measurement, allowing you to document additional information, such as dietary changes or other relevant observations.

Data on the child are presented on percentile grids in relation to age and body weight. This solution allows for ongoing monitoring of the growth rate and comparison of results with population norms.

The application is informative and auxiliary. The results obtained do not constitute a medical diagnosis or medical advice. If you have any doubts about your health condition, you should consult a qualified medical professional.

07. CLEANING AND STORAGE

Clean the appliance with a soft, slightly damp cloth. When cleaning is complete, wipe the surface dry thoroughly.

Avoid excessive contact of the device with moisture, as it may enter the interior of the enclosure and cause damage to electronic components.

Do not use chemicals such as emulsions, lotions, cleaning pastes or other preparations containing aggressive substances for cleaning. Do not use abrasives that may damage the surface of the appliance.

When not in use, store the appliance in a cool, dry place, out of the reach of children.

Attention! Failure to follow the use, cleaning, and storage instructions set forth in this manual may result in damage to the device, malfunction, and a shortened expected service life.

08. SPECIFICATION

Minimum Weight Measurement: 0.2kg

Maximum load: with pad: 20 kg/without pad: 180 kg

Measurement accuracy: 10 g

Units of measurement: kg/lb
Weight Color: White
Overlay Color: Multicolor
Display: LCD 73 x 30 mm
Device dimensions: with cap: 53 x 27 x 6 cm/without cap: 30 x 30 x 1.5 cm
Material: Tempered glass, ABS
Power supply: 2 x CR2032 (3 V DC) - not included
Overload indicator: Yes
Low battery indicator: Yes
Automatic shut-offs: yes, after 1 min
Display rotation: Yes
- Tare option: yes
- HOLD Option: Yes
Device weight: with cap: 1.93 kg/without cap: 1.3 kg
Operating Temperature Range: 5-35°C
Package dimensions: 34 x 13 x 34.7 cm
Package weight: 2.23 kg
BT connectivity: yes
Dedicated App: Neno
Presentation of analytical data: access to full parameter analysis in the application
Newborn data presentation: body weight, percentile grids
Frequency Range: 2402–2480 MHz
Maximum radio frequency power emitted in the frequency range: 2.84 dBm

09. WARRANTY CARD

The product is covered by a 24-month warranty. The warranty conditions can be found at:

<https://neno.pl/gwarancja>

Details, contact and service address can be found at: <https://neno.pl/kontakt>

Kit specifications and contents are subject to change without notice. We apologize for any inconvenience.

KGK Trend declares that the Neno Integro device complies with the essential requirements of Directive 2014/53/EU.

The text of the declaration can be found at the link: <https://neno.pl/download/DOC/deklaracja-CE-Neno-Integro.pdf>

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